

WATER COOLING SYSTEM

Bringing water cooling to the Dynamic and Ecomak ranges.

- Fine spray mist directed onto the cutting surface, reducing heat and dust.
- Helps to prevent thermal damage to the dental tissues.
- Especially useful when widening diastemata or making reductions to large overgrowths.
- Clever design allows the fine mist to be expelled in the correct place even on the curved handle when the head has been positioned in a different orientation.



Water cooling handpieces are compatible with our whole range of heads and burs



Head rotates into four different positions




SCAN ME

For introductory offers and more information on our water cooling system, please scan this QR code for more details

Water cooling mist spray handpiece

Easy to replace water pipe using simple connections

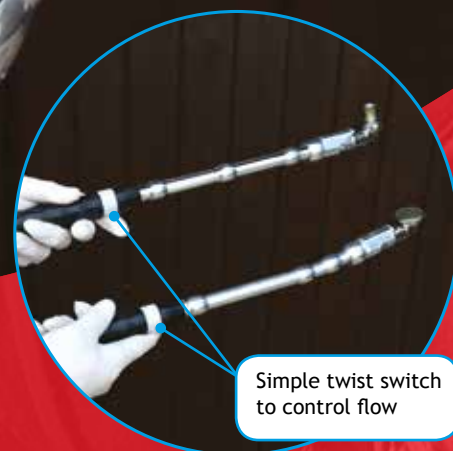

Simple twist switch to control flow

A revolution in equine dentistry

- ✦ Developed to be the most user-friendly dental system on the market
- ✦ Dynamic adjustable torque motor - new & improved
- ✦ Cross/Key/Square/Slot Drive Cable options
- ✦ Lightweight aluminium handpieces
- ✦ Longer lasting heads
- ✦ Compatible with other power systems
- ✦ Reduced servicing costs & option for DIY servicing
- ✦ Available in a variety of complete kits to suit the different needs of practitioners
- ✦ Comes in motor on the belt or hanging motor version
- ✦ Made in the UK

DYNAMIC Motor Powerful, quiet, brushless motor. Adjustable torque system. The motor unit contains the on/off switch and the variable speed dial, torque control and foot pedal dock.

DYNAMIC Water Cooling Kit
All Dynamic kits are available with water cooling upgrade (see front cover of this brochure for details).

NEW



Carry case

DYNAMIC Cross Drive
Unique DYNAMIC Cross Drive. Cuts vibration for a smoother power delivery. Also available in our own Key Drive, Square Drive, Slot Drive and Shear Pin Drive cable head fittings.

Long straight handpiece with DIVERSE head and Applecore burr

Long curved handpiece with Disc head

Battery charger unit



Grease gun GREASEGUN
Dynamic food grade grease GREASE

Disc/Burr Removal Kit




DYNAMIC Base Kit (Hanging Motor and Belt Motor Option)

Includes belt, or bracket, Dynamic motor, one Makita 40 volt 2.5 ah battery, Makita battery charger. (Braces and Holster are optional extras).
DKITBASE
DKITBASEH

DYNAMIC Battery:

The Dynamic belt holds one 40 volt Makita Battery. The supplied 2.5 ah battery (treats up to 10-12 horses) is lightweight very durable and very powerful. 4 ah batteries are available for busy practitioners (treats up to 18-20 horses).

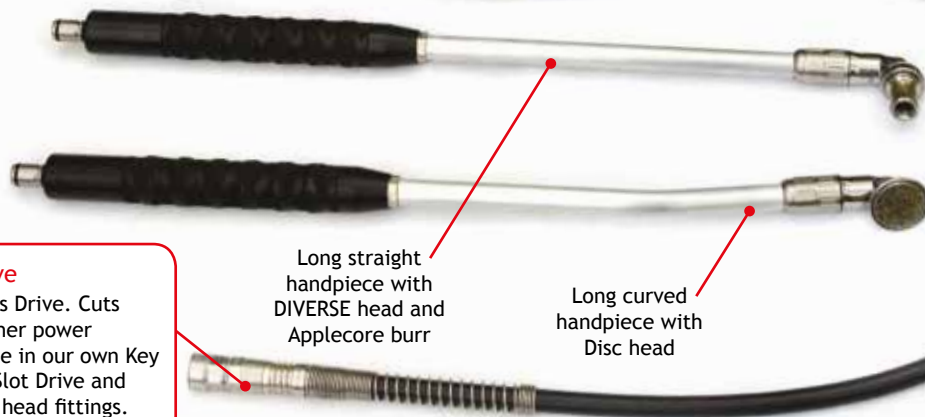

DYNAMIC Complete Kits

Motor on Belt and Hanging Motor

DYNAMIC
DENTAL EQUIPMENT

3

DYNAMIC Disc and Diverse Kit

Great for busy practitioners
As DYNAMIC Disc Kit plus:
Straight Dynamic Diverse
handpiece with
a medium applecore burr and
spanner.
Straight shaft option also
available.
DKITCOMP



WATER COOLING UPGRADE
AVAILABLE (SEE FRONT COVER)

DYNAMIC Disc Kit

Great for everyday dentistry

This kit includes belt, Dynamic motor,
Makita charger, one 40 volt 2.5 ah battery
and the Dynamic drive cable.
Curved Dynamic Disc handpiece. Also
included Disc Removal Kit, grease gun,
Dynamic food grade grease and basic
Dynamic carry case.
Straight shaft option also available.
DKITDISC



WATER COOLING UPGRADE
AVAILABLE (SEE FRONT COVER)

Hanging Motor Kits

Includes bracket, Dynamic
motor, one Makita 40 volt 2.5 ah
battery, Makita battery charger.

DKITDISC
DKITDUOH
DKITCOMPH

Also available in Disc,
Duo and Disc & Diverse
options (as per the motor
on belt kits)



Foot Pedal
recommended
as an optional
extra.



WATER COOLING UPGRADE
AVAILABLE (SEE FRONT COVER)

DYNAMIC Duo Kit

Save money, both heads, one shaft
As DYNAMIC Disc Kit plus:
Dynamic Diverse head (no shaft) with
a medium applecore burr and spanner.
Straight shaft option also available.
DKITDUO



WATER COOLING UPGRADE
AVAILABLE (SEE FRONT COVER)

DYNAMIC Accessories (Optional Extras)

**DYNAMIC
Basic Carry
Case**
DCASETB



**DYNAMIC Battery Charger and
Spare Batteries** are available to
purchase separately. The battery
is available in different Ah
depending on the requirements
of the practitioner.
MBATCHS40
MBAT2.5-40V



**DYNAMIC
Trolley
Case**
DCASE



DYNAMIC Holster
Ergonomic way to
carry your handpiece
and drive cable.
DHOL



DYNAMIC Foot Pedal
This is a plug in option ideal for
more intricate procedures. Comes
with cover. DFOOTMAK

DYNAMIC Braces
Hard wearing, wipe clean material
designed with maximum comfort in
mind. DBRACE



Unique. Cuts vibration for a smoother power delivery.



Heavy Duty Cable

Best seller

DDRIVE120

Cross-drive



DRIVE CABLES

Available in 120 cm and 160cm lengths with a 6mm reinforced inner cable. Available in cross, slot, key or square drive. Please contact us or your distributor to discuss your specific requirements.

LONGER LENGTH AVAILABLE

DDRIVE120P
DDRIVE160P



Heavy Duty PLUS Cable
Our most robust cable for heavy users.

Cross-drive

NEW ROBUST CABLE



Slot-drive

EDRIVE120
EDRIVE120P
EDRIVE160P



Key-drive

KDDRIVE120
KDDRIVE120P
KDDRIVE160P



Square-drive

SDDRIVE120
SDDRIVE120P
SDDRIVE160P

Sectioning Handpieces, Head and Burr

- 90 degree head.
- Used to section teeth during difficult extractions, reducing the need for more invasive extraction techniques.
- Low profile for maximum accessibility and manoeuvrability.
- Head compatible with all dynamic shafts and handles.
- Recommended to be used with water cooled handpieces.
- Used in conjunction with carbide sectioning burs.
- Available in cross drive, key drive, slot drive or square drive.



NEW

For more details and prices



DENTAL DRILL 75MM LENGTH 3.175MM DIA
DENTAL DRILL 50MM LENGTH 3.175MM DIA



NEW

DYNAMIC Water Cooling Handpiece
All Dynamic kits are available with water cooling upgrade (see front cover of this brochure for details).



Handpiece, Head and Burr Combinations



Dynamic disc short grip handpiece with disc head
DDISCSG

Dynamic short grip handpiece with drumhead
DDRUMSG


Dynamic incisor disc short grip handpiece
DINCISORSG

Dynamic incisor cut off disc short grip handpiece
DMANDRELSG

Dynamic Disc straight medium handpiece
DDISCS(MEDIUM)



Dynamic Disc curved long handpiece
DDISCC(LONG)

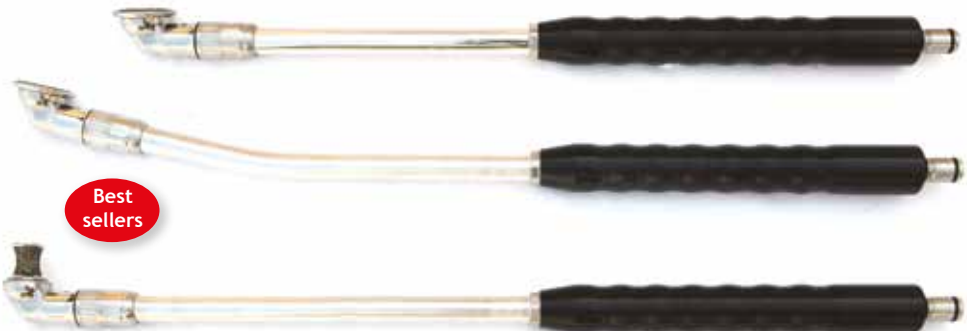

Best sellers

Dynamic Diverse straight long handpiece D90S(LONG) with medium applecore burr
DAPPLEM

Dynamic Diverse straight medium handpiece with buccal burr
D90S(MEDIUM) DBURRB

Dynamic Drum straight medium handpiece
DDRUM(MEDIUM)

Dynamic Diverse curved long handpiece with diastema burr
D90C(LONG) DDIASST(SLIM)

All handpiece shafts are available with straight, curved, medium or long shafts

DYNAMIC Diverse Burr options

Can be used with any Diverse handpiece and will allow for a number of treatments including buffing canines, profiling of the 06's and treating diastemata.

Disc/Burr removal kit -



Short/Extra Slim Diastema Burr 20 x 3 mm
DDIASTXS



Short/Slim Diastema Burr 20 x 4 mm
DDIASTSLIM



Short/Wide Diastema Burr 20 x 5 mm
DDIASTWIDE



Medium/Slim Diastema Burr 25 x 4 mm
DDIASTSLIM-25X4



Long/Slim Guarded Diastema Burr 30 x 4 mm
DDIASTSLIM-30X4



Cone Burr 14mm x 5mm
DBURRC



Small Applecore Burr 14mm
DAPPLES



Medium Applecore Burr 18mm
DAPPLEM



Medium Applecore Burr with Canine End 18mm
DAPPLEMC



Large Applecore Burr 22mm
DAPPLEL



Medium Squat Burr with Canine End 18mm
DBURRMC



Buccal Burr 25mm x 25mm
DBURRB

The 'New' ECOMAK system is a cost effective solution developed for starter/occasional practitioners, or those looking for a spare back up motor

- Key/Square/Slot/Cross Drive Cable options
- Supports the use of the revolutionary DYNAMIC Handpieces
- Lightweight
- DIY Servicing
- Made in the UK



Carry case
ECASET B

ECOMAK Motor

Powerful, brushed, battery driven. The motor has an on/off switch and speed dial at the top of the housing. It also has a hanger loop at the base of the motor in case the practitioner prefers to use without a belt.

Battery

Long straight handpiece with DIVERSE head and Applecore burr

Long curved handpiece with Disc head

DYNAMIC Cross Drive
Unique DYNAMIC Cross Drive. Cuts vibration for a smoother power delivery

DYNAMIC Water Cooling Kit
All ECOMAK kits are available with water cooling upgrade (see front cover of this brochure for details)

NEW



MKITBASE

ECOMAK Belt MBELT

Heavy duty, yet lightweight. The spring loaded motor connector rotates allowing the motor to move in the direction of the flexible cable.

ECOMAK Holster (Optional)

Grease gun, Dynamic food grade grease
GREASEGUN
GREASE

Disc/Burr Removal Kit

Battery charger

ECOMAK Battery

18 Volt Lithium Ion 5Ah Makita battery which can treat up to 4 horses per battery. The battery takes 45 minutes to charge. Extra batteries can be bought from any DYNAMIC/ECOMAK distributor or from any Makita Stockist, helping to reduce problems when travelling or shipping abroad.



ECOMAK Kits

ECOMAK Disc Kit



Belt, ECOMAK motor, 18 volt 5Ah Makita Battery and a single unit battery charger, Dynamic drive cable, curved Dynamic Disc handpiece. Also included Basic repair tools, grease gun, Dynamic food grade grease and basic ECOMAK carry case.

Straight shaft disc option available
MKITDISC



All ECOMAK kits are available with water cooling upgrade (see front cover of this brochure for details)

ECOMAK Duo Kit



Disc Kit plus: Dynamic Diverse Head (no shaft) with a medium applecore burr and spanner.

Straight Shaft option available.
MKITDUO



All ECOMAK kits are available with water cooling upgrade (see front cover of this brochure for details)

ECOMAK Disc and Diverse kit



Disc Kit plus: Dynamic Diverse handpiece with a medium applecore burr and spanner.

Straight shaft disc option available

MKITCOMP



All ECOMAK kits are available with water cooling upgrade (see front cover of this brochure for details)

Optional Extras



ECOMAK Holster

Strong and simple, ideal for storing handpiece mid-way through dental procedures.

MHOL

ECOMAK Braces

Hard wearing material, wipe clean. Designed with maximum comfort and breathability in mind.

MBRACE



ECOMAK Battery Charger and Spare

Batteries are available to purchase separately. The battery is available in different Ah depending on the requirements of the practitioner.

ECOMAK Foot Pedal

This is a plug in option ideal for more intricate procedures. Comes with cover.

DFOOTMAK



Servicing and Aftercare

Servicing your DYNAMIC or ECOMAK powered equipment is an essential part of the after-care required to ensure that it continues to function effectively and efficiently.

The service options available:

- **DIY Servicing**

Use the DYNAMIC online guides available on the website for reference, with the DYNAMIC Disc/Burr removal kit (available to buy from your local distributor).

- **DYNAMIC Equipment Service Centre**

It is recommended that all complete dental power tool systems (DPTS) receive an annual service in order to prolong the life span of the equipment. Regular servicing should also reduce long-term servicing costs.

For a very busy practitioner treating 1000 - 2000 equines a year, it is recommended that DPTS is serviced twice a year.

A number of DYNAMIC global distributors have trained in-house technicians to service the DYNAMIC and ECOMAK powered equipment, and will turnaround this equipment as quickly as possible and will always, where possible, discuss any work that needs to be undertaken, and the costs involved prior to starting the service. The DYNAMIC Equipment Service Centres are also able to work on any repairs to other brands of equipment

To find your local Service Centre contact us on

0044 (0) 1934 287818

info@dynamicdentalequipment.com

www.dynamicdentalequipment.com/servicing

Speculums and Accessories

Hard wearing speculums, made to our in-house design

Mini Super Deluxe Pony Speculum

SPEC-PONY
Not Available in USA



This Mini SD speculum has been designed to create a better overall experience for our small equine patients. It has smaller opening increments and bite plates which makes it more comfortable. The bite plate design also creates a slight spreader effect making it easier to work.

Offset Gunther Speculum

SPEC-GUNTEROFF

Screw handle opening ratchet with 9 offset positions available to assist with difficult oral extraction. Works well with Advanced MTE kit and our MIR kit

NEW

Comes with Biothene Straps and Rubber Bite Plates. Easy to clean.



Rubber



Extra Wide Bite Plates with rubber lining to prevent chipping. Replacement inserts sold separately
BITEPLATERUBBER



Extra Wide Bite Plates with spreader effect
BITEPLATEXW



Parrot Plates to fit equines with extreme incisor congenital abnormalities. Can be sold individually or as a pair BITEPLATERUBBERP



Padded Gum Bars
Stainless Steel, extended Gum Bars
BITEPLATEGUM

Rubber

Sold individually, either Flat or Toed options available. Screws sold separately.



Replacement Inserts
BITEPLATERUBINSERT



Cheek Retractor
CHEEKRETRACTOR



Biothene Strap Set
SPEC-STRAPS



Incisor Speculum
SPEC-INCISOR

X-Ray Speculum

SPEC-XRAY

Hands free. Offsets the lower jaw to the left or the right. Gives an open mouth view. Can be used to investigate roots, view the occlusal surface of the cheek teeth or get offset view of the dorsal-ventral

Offset position



X-Ray Bite Plates

BITEPLATEXRAY

Specialist bite plates for equine head radiographs. Only intended for taking X-rays of incisors. Fits most speculums on the market.

Dynamic Headlight with USB/Rechargeable Battery

LIGHT-HEAD-USB

BEST SELLER

The New Improved Dynamic headlight is incredibly bright with 6700k superclear daylight 1300 Lumen XML2 LED. It has a long life, waterproof rechargeable USB Battery. Incredibly lightweight but robust with 27mm reduced diameter stainless steel housing and a comfortable heavy duty headband.




Speculums and Accessories

Hard wearing speculums, made to our in-house design

DYNAMIC
DENTAL EQUIPMENT

9

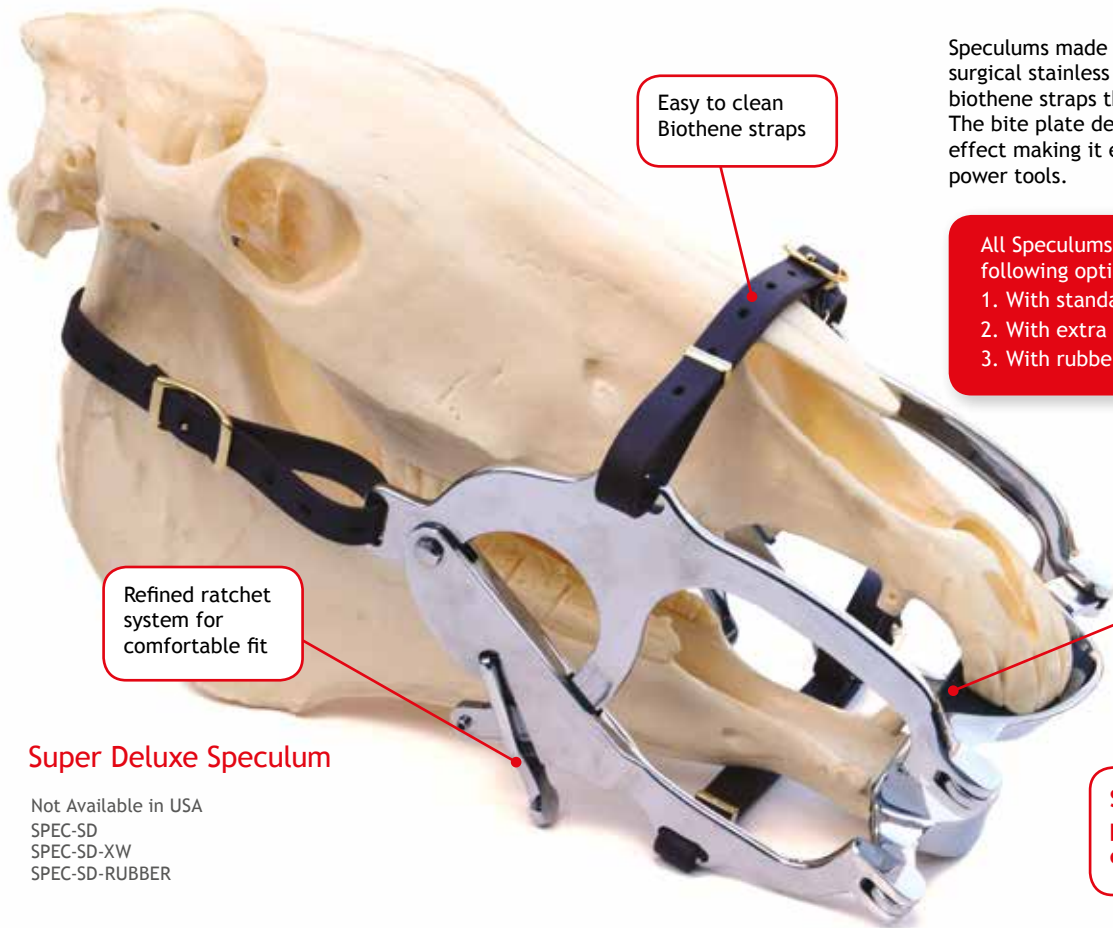

Easy to clean
Biothane straps

Speculums made from top quality polished surgical stainless steel with hardwearing biothane straps that are easy to clean and use. The bite plate design creates a slight spreader effect making it easy to work, especially with power tools.

All Speculums are available with the following options:

1. With standard bite plates
2. With extra wide bite plates
3. With rubber lined extra wide bite plates

Refined ratchet
system for
comfortable fit

Bite plates
angled to
accommodate
incisor angle

Super Deluxe Speculum

Not Available in USA
SPEC-SD
SPEC-SD-XW
SPEC-SD-RUBBER

Servicing and spare
parts available. Please
contact us for details

SE Speculum

Spreader effect
with five ratchet
settings for a
good fit.



SPEC-SE
SPEC-SE-XW
SPEC-SE-RUBBER



Millennium Speculum




SPEC-MILLENNIUM
SPEC-MILLENNIUM-XW
SPEC-MILLENNIUM-RUBBER

Spreader effect and 15 small ratchet settings for an exact fit.



Speculum Spreader Bars available

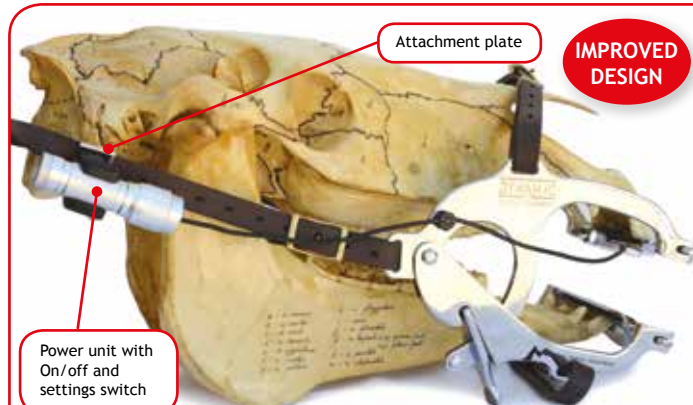

Attachment plate


**IMPROVED
DESIGN**

Power unit with
On/off and
settings switch

DYNAMIC Magnetic Speculum Light

SPEC-LIGHT-DYNAMIC

A powerful and durable magnetic light that fits onto speculum bite plates using ultra strong magnets. Super bright clear daylight LED with 3 brightness settings. Comes with rechargeable water-resistant long-lasting battery and charger.

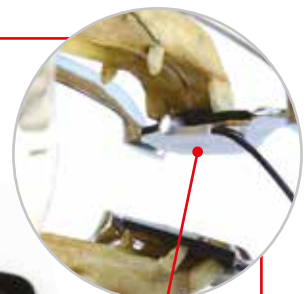

Powerful LED
with super slim
and ultra strong
magnets



DYNAMIC Full HD Intra Oral Camera ENDOFULLHD

Great quality at a very competitive price.

- 5.0 mega pixel macro lens auto focus (1.5cm - 3.0cm) camera module image resolution
- 10 bright LEDs
- Water resistant strong easy clean lightweight alloy body (600mm length)
- Image storage to PC or Android phone or tablet
- Hard case included
- 90 degree wide viewing angle
- Image flip / reverse

NEW

Robust and fit for purpose

Super strong carry case (not to scale)

Long lead to connect to PC, with software supplied on USB stick

Incredible high definition images

Easy to use and set up

Water-Resistant

*PC not included

Great quality at an affordable price

DYNAMIC are proud to be the PowerFloat European Distributor and Service Centre.

The PowerFloat is available in both Battery or Mains Power form.

We stock a wide selection of 'Quick Connect' attachments, burs and discs.

For more information about the complete range and current prices please contact us.

PowerFloat



NEW

Each variety of PowerFloat can be bought on its own or as part of a Basic or Complete kit.

A **Basic Kit** will include:

- Pelican® Carry Case (Lifetime warranty)
- PowerFloat Wheel Removal tool
- LED Headlight
- Grease Gun
- Grease Tube: 15oz tube
- Two Nylon Mouth Gags (large and small size included)
- The Guarded Right Angle shaft with a diamond disc

A **Complete Kit** includes all the same as the Basic Kit but has the addition of:

- Diamond Grinding Wheel
- Long Rotary Attachment
- Choice of Short Guarded Right Angle (w/ Diastema Burr) or Short Un-Guarded Right Angle (w/ Apple Core Burr) Attachment
- Un-Guarded Right-Angle Attachment (w/ Chamfer Burr)
- Round Nose Burr (3/8")
- Upper 11s Hook Burr
- Lower 11s Hook Burr



For more information and prices

Minimal Invasive Buccotomy, Transbuccal Screw Extraction

The MTE kit allows a safe technique for the removal of roots, fractured teeth or those that have insufficient clinical crown length. We recommend the following equipment should be used in conjunction with this MTE kit:

- 1) The Offset Gunther Speculum:- enables better access to the cheek teeth.
- 2) Duckbilled Forceps for use during the intraoral screw extraction technique.

Made in Germany

Standard Kit MTE-STOLL-STANDARD

Contents: Tap M6, Trocar sharp long, Trocar sharp short, Trocar round long, Trocar round short, Drill sleeve short serrated, Drill sleeve long serrated, Fragment forcep with spring handle, Elevator for root tips slim, Elevator sharp spoon, Elevator symmetric, Elevator asymmetric, Threaded pin (screw) M6, Autoclave sterilisation container with silicone mat, Metal slotted hammer, Drill 3,2 mm, Drill 5 mm, x 200mm, Intraoral screws 6mm

Advanced Kit MTE-STOLL-ADVANCED

Contents: Same as Standard kit plus Tap M4, Tap M7, Threaded pin (screw) M4, Threaded pin (screw) M7



All components are available seperately - see website for details

Absorbent Paper Points

Sold as a variety pack and as single sizes. There are a range of lengths and widths to suit your Equine Endodontic treatment.

Also available in the Pferdfit-Dental range

Low and High Speed Burrs

A variety of lengths and diameters of ball shaped, Lindem burrs and Cindem cylindrical burrs. Low speed burrs to fit contra angle handpiece, latch type, and high speed burrs for turbine and contra angle handpiece with chuck.

Autoclave Tin

Tin to store Hedstrom files, fully autoclavable. Comes in small and large sizes.

Hedstrom Files

Sold in pairs, highly flexible Hedstrom files in a range of lengths and diameters. Most files are available with Latch type connection to an Endodontic handpiece.



See website for prices


Forceps and Fulcrums

Tools designed to aid successful molar extraction

Heavy Duty Duck Bill Forceps

FORCEPS-DUCKBILL

For use when extracting cheek teeth orally and access within the mouth is limited. The duck bill is designed to enable the practitioner to get further purchase, compared to other heavy duty forceps.

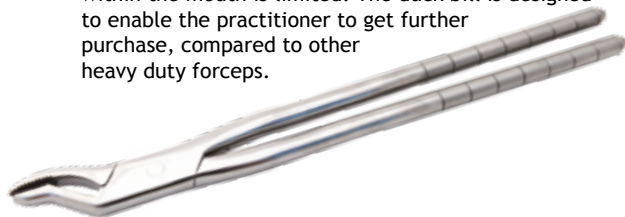

Root Fragment Forceps

FORCEPROOT

Stainless Steel 42cm Root Fragment Forcep with 60mm serrated head for removing difficult fragments.



Slab Fracture Forceps

FORCEPSLAB

Stainless Steel 42cm Slab Fracture Forcep with 30mm serrated head for removing difficult slab fractures.



Designed in conjunction with Professor Paddy Dixon MVB PhD MRCVS

Long Nose Forcep

FORCEPLONGNOSE 35/50/70

The narrow profile is ideal for root fragment extractions. Length: 43 cms. Three nose lengths available.



Fulcrums



Nylon 6 Tubular Fulcrum Kit

FULCRUM6HEADS

Fulcrum with 6 white tubular nylon blocks, designed to fit snugly on the tooth



Nylon 8 Block Fulcrum Kit

FULCRUM8HEADS

Fulcrum with 8 numbered white nylon blocks. Comes in a sturdy case

NEW



Mandibular Fulcrums

FULCRUM-SET5-MAND

Set of 5 stainless steel fulcrums ideal for cheek teeth extractions (2mm, 3mm, 4.5mm, 6mm, 8mm). Newly designed to eliminate the need to change blocks.

Rostral Caudal Fragment Forcep

FORCEP-ROSTAL

46cm. Tightens to the tooth. Excellent for difficult molar fractures.



Grips fragment front and back

Hewson Style Forcep

FORCEPS15-HEWSON

Ideal for premolar cap removal.

Length: 40 cms



4 Root Forcep

FORCEPS16-FOURROOT



Pony Forceps

Stainless steel forceps designed to fit the smallest equine teeth.

New lower profile spreader



Spreader




Serrated



FourRoot

FORCEP-PONY-SPREADER/ SERRATE/4ROOT

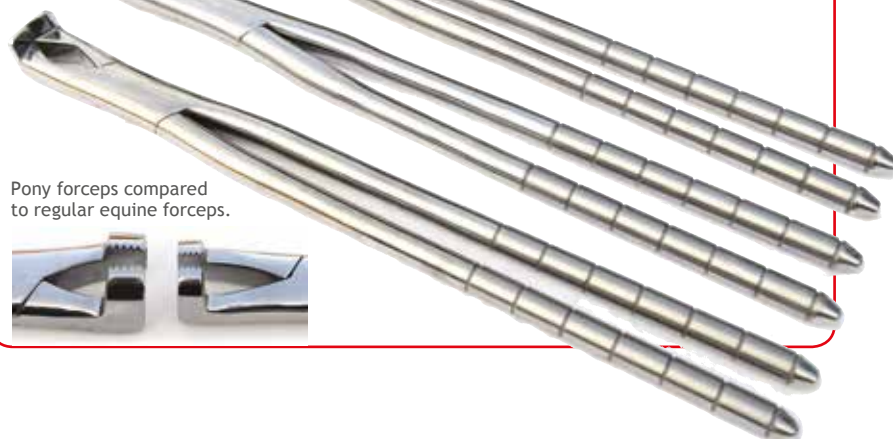

Pony forceps compared to regular equine forceps.



Multi-Pick T Handles and Forceps

Fully interchangeable to create the perfect pick for any job

DYNAMIC
DENTAL EQUIPMENT

13

A range of fully interchangeable handles and inserts that can be purchased individually or as a set. Each Insert Head slots easily into any of the Multi-Pick Handles and tightened into any rotation with a 10mm spanner and allen key.

Starter Kit MULTI-PICK-SET1

Universal Handle, 0° Handle, Spanner, Allen Key, Screw, Zip-up case and 6 insert heads: Small (Alveolar, Elevator, Spreader), Medium (Elevator, Spreader) and a Large Spreader.



Complete Kit MULTI-PICK-SET2

Full range of 4 Handles and 13 insert heads as well as Universal Handle, Spanner, Allen Key, Screws and Zip-up case



Stronger Inserts

Stronger Design



Alveolar
MULTIA - S/M/L



Spreader
MULTIS - S/M/L



Curette
MULTIC - S/M/L



Elevator
MULTIEL - XS/S/M/L

Insert Heads

Multi-Pick Handle

MULTI-PICK-THANDLE0 & MULTI-PICKTHANDLE 7/15/25

Four handles each with an angled head 0°, 7°, 15°, 25°



T Handle Excavators - Set of 2 (35mm head length)

EXCAVATORS-SET2

Left and Right. T-handle for greater leverage and control.

A zip up case is included.

NEW



DYNAMIC "Pro" Forceps

Introducing a new and improved range of DYNAMIC "Pro" Forceps which have been designed in collaboration with leading equine veterinary dental specialists.

These forceps have been developed on a heavy duty frame but are lightweight enough to allow for one handed operation to ensure comfort and reduce operator fatigue.

This new design significantly increases the success rate of oral extraction. They have a square box joint which gives the forceps added strength and also helps prevent them from slipping off fulcrums during extractions.

Length: 48cm

Universal Spreader Forceps, available in 3 sizes as well as 2 serrated options:



Slim Spreader Forceps (A) 4mm
FORCEPS-PRO-SPREADER-S

Thinner/wider wedge head design to aid initial spread



Medium Spreader Forceps (B) 5mm
FORCEPS-PRO-SPREADER-M



Wide Spreader Forceps (C) 6mm
FORCEPS-PRO-SPREADER-W



3 Root Inside Forceps (G)
FORCEPS-PRO-3ON -



4 Root Forceps (F)
FORCEPS-PRO-4ROOT



3 Root Offside Forceps (H)
FORCEPS-PRO-3OFF



Concave (I)
FORCEPS-PRO-CONCAVE



Universal Serrated Forceps (D)
FORCEPS-PRO-SERRATED-U
Suitable for both Mandibular and Maxillary cheek teeth



Maxillary Serrated Forceps (E)
FORCEPS-PRO-SERRATED-M
Wider jaws for use on the wider upper maxillary cheek teeth, allow the handles to be closer together for comfortable one handed operation.

Multi-Pick T Handles and Forceps

Dentalvator

Revolutionary powered equine dental elevator.

- Different sized replacement elevators available.
- Reciprocating blade to elevate the periodontal tissue.
- Easy to remove and reposition the blades.
- Suitable for incisor, canine, wolf teeth and cheek teeth.
- Adjustable grip handle for both incisor and cheek tooth work.
- Available in cross, key, slot or square drive.
- Ground breaking new extraction technique.
- Much more comprehensive elevation leading to many more successful extractions.
- Greatly reducing time spent elevating prior to extraction.

NEW

Dentalvator Kit

Includes: Molarblades set (1 Small, 1 Medium, 1 Large Blade)
VATOR

Pack of Angled Roof Blades

VATORB-A

Wolf / Canine Incisor Blade

VATORB-W

Pack of Molarblades

VATORB-M

MIR Kit MIRKIT/MIRKITL

For minimally invasive repulsion of dental fragments or fractured teeth. For use by vet surgeons in conjunction with careful radiographic positioning through a small trephine or drill hole. Six pins. Each has a block end that slots into the hammer to aid removal should the pin become lodged in an incorrect position.
(Replacement pins available)

Developed in conjunction with Chris Pearce MRCVS from Equine Dental Clinic



New Improved Design



Available in short or long length pins

Universal Dental Hammer

MIRHAMMER


Included in the MIR kits, but also available separately. Works well with set of 6 elevators and set of 4 short grip elevators. Used to aid wolf/canine/incisor extractions.



Extended Endo File Holder FILEHOLDER

40cm handle with knurling for easy grip. Head can be adjusted to suitable angle. Will hold universal sizes of H and K files (Not included)


New Improved Design

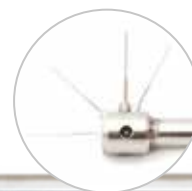

Universal K files pack of 5
FILEK-PACKOFFIVE



Diastema Putty

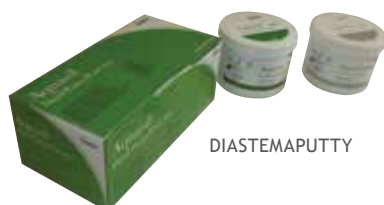

DIATEMAPUTTY

Fast setting impression material for diastema packing. 450 ml catalyst and 450ml base.

Extended General Nerve Block Syringe

40cm from base to needle

Ideal for administering intra-oral local anaesthetic for mandibular/general nerve blocking prior to molar extraction.

Both available with 3 different syringe holder ends:

- 1.8ml UK Thread
- 1.8ml Euro Thread
- 2.2ml UK Thread

NEW

SYRINGELONG

Wolf Tooth Nerve Block Syringe

25cm from base to needle

Extended Intra-Oral LA Syringe. Ideal for administering local anaesthetic for nerve block prior to wolf tooth extraction.

SYRINGESHORT



Individual Incisor Forceps

FORCEPS-INCISOR-DUCKBILL

28 cm Duck Bill Incisor Forceps

Ideal design for extracting incisor teeth.

Unique angle allows better access for secure purchase on the tooth.



Set of 6 Incisor Elevators

ELEVATOR-SET-6

Set of 6 incisor elevators with nylon handles and stainless steel base designed to permit gentle hammer usage.

- 3 different angles allows instruments to follow the natural curve of the teeth.

- Comes with 8mm & 12mm head widths to suit different widths of incisor teeth



Set of 4 Short Grip Elevators

ELEVATOR-SET-4

Designed for extracting wolf teeth, but can also be used for canines, incisors and root fragments extractions. Included in the set are two 4 mm elevator (one straight and one concave) and two 5mm elevator (one straight and one concave). The concave curved tip hugs the palatal aspect of the wolf tooth reducing the risk of iatrogenic trauma to the palatine artery.

Cold sterilisation only



Cold sterilisation only

Basic Burgess Set

BURGESS-BASIC-SET

Tough matt handle grip with two stainless steel screw-in extensions (one straight, one offset) and four cut-out elevator heads (6mm, 8mm, 10mm, 12mm).



Sturdy case included

Complete Burgess Set

BURGESSSET-COMplete

Stainless steel 26cm wolf tooth forceps, two stainless steel elevators, tough matt handle grip with two screw-in extensions (one straight, one offset) and four cut-out elevator heads (6mm, 8mm, 10mm, 12mm).

Burgess Heads (available in the following sizes):

Cut-out head sizes 6mm, 8mm, 10mm and 12mm.

Complete head sizes 4mm, 6mm, 8mm, 10mm, 12mm and 14mm



Sturdy case included

Long Elevator Kit

ELEVATOR-OFFSET

Two stainless steel elevators with fine heads, offset for easy use with a speculum.

Comes in a canvas case.

Length: 28.5cm




Compound Rongeurs

COMPOUNDROUGEURS

Stainless steel with compound action and angled head.

Length: 24cm

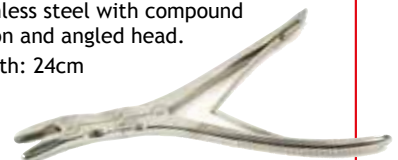

Ideal for wolf teeth and fragment extractions plus deciduous incisor cap removal.

Curved 21cm Wolf/Incisor Tooth Forceps

WOLF8.5-CURVED

Stainless steel 21cm Curved forceps. Great for wolf teeth and incisor cap extractions, fragment removal, as well as removing plaque from the canines.



Wolf/Incisor Tooth Forceps

(Also available in a set of three)

Stainless Steel 18cm forceps with easy grip knurled handles. Great for wolf teeth and incisor cap extractions, fragment removal, as well as removing plaque from lower canines.



WOLF7-SET-3

WOLF7CURVED

WOLF7SPREADER

WOLF7OFFSET

Picks, Probes and Mirrors

Instruments to compliment your kit

Stainless Steel Picks

New design - coloured bands for easy diagnosis

Banded Probe PICK-BANDEDPROBE

Periodontal tissue probe. Measuring depth of perio pockets.

Equine Explorer PICK-EXPLORER

For pulp probing.

Gingiva 90° Pick PICK-GINGIVA

General pick - cleaning diastemas.

Dental Scaler Excavating PICK-DENTALSCALER

Used for excavating occlusal tissue.

Mini Spatula MINISPATULA

NEW Replacement Screw-in Tips

Diagnostic Kit DIAGNOSTICKIT4

Mirror, Alligator Forceps, Explorer Probe, Dental Scaler, Banded Periodontal Probe, Gingival Pick and zipped case.



Fantastic
saving

Peridontal Forceps

FORCEPS-PERIO-30

FORCEPS-PERIO-35

High quality 30cm long stainless steel forceps with very fine tips (also available in 35cms). Designed to remove forage that accumulates in periodontal pockets and valve diastema between equine cheek teeth)

Long Alligator Forceps

FORCEPS-ALLIGATOR-STRAIGHT

Right angled head for removal of food and Hay particles trapped between teeth.

Length: 35cm

we also have a 30cm version

New
Stronger
Design

Straight
handle

Mirrors

Dental Mirror with 40mm lens (41cm) MIRRORM

Best
seller

Magnifying Dental Mirror with 40mm lens (41cm) MIRRORMAG

Small Mirror with 25mm lens (41cm) MIRROR25

Replacement mirror glass available



DIAMONDFLATFILE

Diamond Files

Stainless Steel with diamond grit. To buff sharp canines and bitseats as well as finishing off incisor realignments. Flat or S-shape file option

DIAMONDSFLOAT



Headstands and Irrigation

Stable, adjustable and easy to assemble

DYNAMIC
DENTAL EQUIPMENT

17

Economy Headstands

Additional tops available separately



Padded Cradle
for everyday
dentistry

HEADSTAND-
BUDGET2

Large,
weighty base
for maximum
stability



**Padded
Rectangular
Rest**
Ideal for
sedated
equines

HEADSTAND-
BUDGET

Pony Headstand

HEADSTANDPONY

Column available
separately
HEADSTAND-
BUDGETC



Hanging Halter HALTER

Padded halter hangs from stocks
and rotates ergonomically
to support sedated equines.
Compact and ideal for travel.



Headrest Cover

HEAD-COVER-ECONOMY OR
RECTANGULAR
Easy clean, removeable cover
for headrest.
Rectangular or cradle fit.

*Fits most
headstands*

Dose Syringe

400ml capacity
DOESYRINGE



*Rounded nozzle
for safe flushing*

Drench Gun

250ml self filling
DRENCHGUN



Silicon Lubricant 30ml available separately

Stainless Steel Wand with robust handle and
fine tip nozzle.
WANDHOZELOCK



*Fits a Hozelock
connector*



Irrigation Pump

IRRIGATIONPUMP

Pressurised water pick with
fine tip nozzle for diastema
and periodontal pocket
cleaning.



**Best
seller**



DYNAMIC Aqua Flusher

DIASTEMAFLUSER

Multi functional flusher
designed for everyday oral
flushing, as well as cleaning
out diastema impactions.



**NEW
Improved
design**

Lightweight, compact battery flusher
comes supplied with two different
diastema attachments and an oral
flushing attachment that attach to the
on/off trigger handle.

Super Concentrated

Fast acting, safe and effective for equine dentistry.
Kills 99.99% of bacteria, fungi and viruses. Makes 200
litres and lasts twice as long as other brands.



WASH-MOUTH
WASH-INSTRUMENT

Instrument
Cleaner and
Disinfectant
5 litres
Peppermint
Mouthwash
5 litres

**Protects
against
STRANGLES
AND
RINGWORM**

UK ONLY

Float Kits, Packages and Buckets

Starting up? Save money with our kits and packages

Premium Float Starter Kit FLOAT-KIT-STARTER3



3 essential premium float handles fitted with 3 tungsten carbide blades and a sturdy roll-up carry case.

Case also available separately

Premium Float Advanced Kit FLOAT-KIT-COMPLETE



6 premium float handles and 6 tungsten carbide blades, sturdy roll-up carry case.

Start up Basic Package



STARTUP1-BASIC

SE Speculum, 3 premium range float handles fitted with 3 tungsten carbide 12, 14 & 16mm, blades.

Start up Plus Package



STARTUP2-PLUS

SE Speculum, 3 premium range float handles fitted with 3 tungsten carbide 12, 14 & 16mm, blades. Plus a Dental Mirror, a Drench Gun and a Diamond S Float.



Economy Bucket with Brush or without Brush
 BUCKET-ECONOMY-BRUSH
 BUCKET-ECONOMY

Start up Complete Package



STARTUP3-COMPLETE

SE Speculum, 3 premium range float handles fitted with 3 tungsten carbide blades. An economy Headstand and a bucket with a brush fitted. Plus a Dental Mirror, a Drench Gun and a Diamond S Float.



Premium Bucket with Brush or without Brush
 BUCKET-PREMIUM-BRUSH
 BUCKET-PREMIUM



Blade Savers

Plastic thick and thin style or stainless steel
 BLADESAVER-PL-THIN OR FAT
 BLADESAVER-SS

Premium Float Handles

Premium Float Handles with easy clean rubber handle. Clearly numbered for quick recognition. The rounded stainless steel shaft is milled flat on one side for correct placement in the mouth. Carefully designed to be free of edges and corners.

Blades sold separately

Regular head use
with bolt-on blade



Slimline head use
with stick-on blade



Incisor 15cm shaft, 5° offset straight head (takes 80mm blade) 10FLOAT		10	-
Premolar 18cm shaft, 30° head (takes 50mm blade) 11FLOAT		11	-
Premolar 23cm shaft 20° head (takes 50mm blade) 12FLOAT		12	22
Premolar/Molar 33cm shaft, 10° head (takes 80mm blade) 13FLOAT		13	23
Lingual and Buccal 38cm shaft, straight head (takes 80mm blade) 14FLOAT		14	24
Upper Back Molar 43cm shaft, 15° upturned head (takes 50mm blade) 16FLOAT		16	26
Premolar Roller Bitseat 23cm shaft, offset head (takes 50mm blade) 17FLOAT		17	-

Premium Pony Floats

Designed for small equines. Stainless steel float handles with rubber handles and low profile heads that take 50mm blades. Clearly numbered for easy recognition.

P-1 17cm 20° Premolar

P-2 23cm Straight, Lingual and Buccal.

P-3 23cm 5° Downturn, Lower Arcade.

P-4 28cm 15° Upturned, Back Molar

Can be sold as a kit (does not include blades)



Premium-Light Float Handles excluding blade

The shorter handle and thinner shaft makes our new range of floats lighter to use. The range still has the great design features of the Premium Range

New Lighter
Weight Design



Pocket Blade 32mm
POCKETBLADESMALL -



Tungsten Carbide rectangular blade
with stainless steel backing plate
or VHB adhesive



Regular finish (80mm or 50mm)



Titanium Nitride Coating
(80mm or 50mm)

See back page for prices

Only the best will do

Cheap blades and inferior regrinds are a false economy...

Rectangular Blades made from sub-micron grade tungsten carbide and cut on CNC machines to produce a superior cutting edge.

Available in **diamond cut** or **master cut** finish.

Best
Sellers

Diamond cut Sharp controlled rasping for everyday dentistry



80mm 80DIAMOND



50mm 50DIAMOND

Master cut Less aggressive for super smooth finishing



80mm 80MASTER



50mm 50MASTER

Diamond cut with Titanium Nitride coating



80mm 80DIAMONDCOATED



50mm 50DIAMONDCOATED

Master cut with Titanium Nitride coating



80mm 80MASTERCOATED



50mm 50MASTERCOATED

Fitted to an aluminium backing plate.
Two size options that fit most brands of float handle.
Available with Titanium Nitride coating to reduce chipping and increase blade lifespan.
Also available without backing plate for slimline float handles

Diamond Cut

Master Cut

Triangular Inserts

Rotate for a fresh razor sharp blade

Triangular inserts (X-files) have three sides and stack together to create three cutting surfaces. The three/four groove pattern is designed to hold a steady line while working on the arcades. These super strong inserts are made from Tungsten Carbide and are designed to be long lasting. The Titanium Nitride coated option increases durability.

Triangular Inserts

Coated Triangular Inserts

X4UNCOATED
X4COATED
X3UNCOATED
X3COATED



Tungsten Carbide Blades

Chip Blades

50mm
80mm



Pocket Blades

32mm
56mm



Oval Blades

58mm



Your local DYNAMIC Dental Equipment distributor is:

ALCYON
BELUX



SCAN ME

For your local distributor location

Images are not to scale. Products may vary slightly from image.

To Find Your Local Distributor: www.dynamicdentalequipment.com