



Surflo®-W

**Intravenous catheter with wings /
and injection port**



**Smooth venous access
Better comfort for the patient**

Material

The catheter is made from radio-opaque, biocompatible FEP. Thanks to the physical properties and the reduced thickness of the catheter the flow rate is optimised.

Internal needle

The quality of the bevel of the internal needle ensures a precise incision.

Conical catheter tip

The perfectly tapered catheter tip allows the catheter to enter the vein without causing any venous trauma and an easy placement of the catheter. In this way the comfort of the patient is guaranteed. The optimal bevel/catheter tip distance reduces the risk of passing through the vein.




Winged catheter hub

The wings ensure easy fixation of the catheter. The EN/ISO coloured hub of the SURFLO®-W with wings confirms the catheter size.

Injection port

This allows the intermittent injection of medicines during infusion. A protective cap using EN/ISO colour coding ensures reliable closure and avoids the risk of contamination. The injection port is latex free.

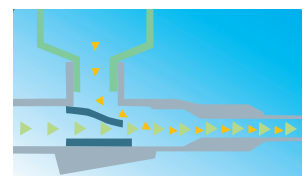

Flash back

The flashback chamber is transparent for easy detection of the flashback. The microporous hydrophobic filter, by allowing the passage of air, ensures reliable closure to avoid any risk of



contamination. The Luer Lock cap allows the extremity of the catheter to be closed after withdrawal of the internal needle.

Obturator

A Luer Lock obturator with the same colour code as the catheter is available for each gauge (except gauge 26). This prevents the formation of clots in the lumen of the catheter due to any interrupted use.




SURFLO®-W with wings and injection port References	EN/ISO colour code	External catheter diameter		Internal catheter diameter mm	Catheter length mm	Flow rate ml/min	Obturator
		Gauge	mm				
SR+DM2619PX	Violet	26	0.64	0.45	19	19	-----
SR+DM2419PX	Yellow	24	0.74	0.55	19	29	SO+DM2419P
SR+DM2225PX	Blue	22	0.90	0.65	25	42	SO+DM2225P
SR+DM2032PX	Pink	20	1.00	0.75	32	59	SO+DM2032P
SR+DM1832PX	Green	18	1.30	0.95	32	103	SO+DM1832P
SR+DM1845PX	Green	18	1.30	0.95	45	96	SO+DM1845P
SR+DM1745PX	White	17	1.50	1.15	45	155	SO+DM1745P
SR+DM1645PX	Gray	16	1.75	1.35	45	225	SO+DM1645P
SR+DM1445PX	Orange	14	2.00	1.55	45	290	SO+DM1445P

SURFLO®-W with wings References	EN/ISO colour code	External catheter diameter		Internal catheter diameter mm	Catheter length mm	Flow rate ml/min	Obturator
		Gauge	mm				
SR+DM2619WX	Violet	26	0.64	0.45	19	19	-----
SR+DM2419WX	Yellow	24	0.74	0.55	19	29	SO+DM2419W
SR+DM2225WX	Blue	22	0.90	0.65	25	42	SO+DM2225W
SR+DM2032WX	Pink	20	1.00	0.75	32	59	SO+DM2032W
SR+DM1832WX	Green	18	1.30	0.95	32	103	SO+DM1832W
SR+DM1845WX	Green	18	1.30	0.95	45	96	SO+DM1845W
SR+DM1745WX	White	17	1.50	1.15	45	155	SO+DM1745W
SR+DM1645WX	Grey	16	1.75	1.35	45	225	SO+DM1645W
SR+DM1445WX	Orange	14	2.00	1.55	45	290	SO+DM1445W



Every day,
as healthcare
professionals,
you are working
to ensure life
and to contribute to
the well-being
of your patients.
Every day,
Terumo Europe
provides you with
the clinical systems,
means and solutions
to keep life flowing.

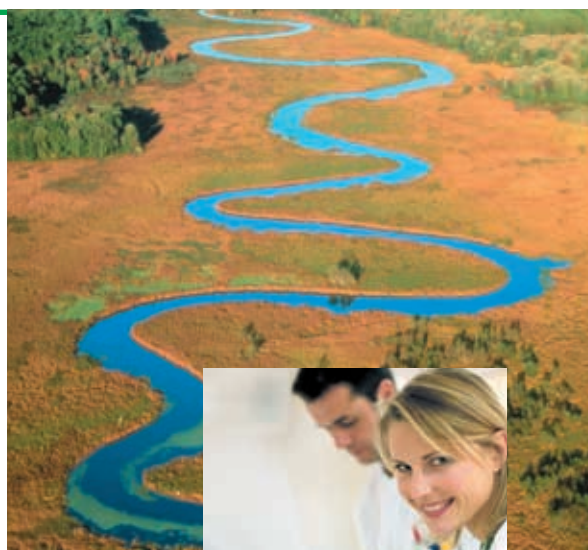

We keep life flowing

TERUMO®

TERUMO CORPORATION
44-1, 2-chome
Hataigaya
Shibuya-ku
Tokyo 151
JAPAN
Tel. +81/3 3374.8111
Fax +81/3 3374.8196

TERUMO EUROPE N.V.
Researchpark Zone 2 Haasrode
Interleuvenlaan 40
B-3001 Leuven
BELGIUM
Tel. +32/16.38.12.11
Fax +32/16.40.02.49

**TERUMO EUROPE N.V.
BENELUX SALES DIVISION**
Researchpark Zone 2 Haasrode
Interleuvenlaan 40
B-3001 Leuven
BELGIUM
Tel. +32/16.38.14.02
Fax +32/16.38.15.55
THE NETHERLANDS
Tel. +0800/023.19.38
Fax +0800/022.34.89

**LABORATOIRES TERUMO
FRANCE S.A.**
1, Parc Ariane
Bâtiment Uranus
Rue Hélène Boucher
78284 Guyancourt Cedex
FRANCE
Tel. +33/1 30.96.13.00
Fax +33/1 30.43.60.85

TERUMO (DEUTSCHLAND) GmbH
Lyonerstrasse 11 a
60528 Frankfurt am Main
GERMANY
Tel. +49/69.66.44.20
Fax +49/69.66.66.826

**TERUMO EUROPE SUCURSAL
EN ESPAÑA**
Avenida de Burgos 16-D
Edificio EUROMOR
Madrid 28036
ESPAÑA
Tel. +34/91 383.90.69
+34/91 383.91.61
Fax +34/1 767.04.60

**TERUMO EUROPE ROME
BRANCH**
Via Simone Martini 143/145
I-00142 Roma
ITALIA
Tel. +39/06.519.61.420
Fax +39/06.503.04.07

TERUMO EUROPE N.V.
Belgien filial Sverige
Nya Varvet
Byggnad 90
SE-426 71 Västra Frölunda
SWEDEN
Tel. +46/31.748.58.80
Fax +46/31.748.58.90

TERUMO U.K.
(Branch of Terumo Europe N.V.)
3 Unity Grove
Off School Lane
Knowsley Industrial Park South
Knowsley
Merseyside
L34 9GT
UNITED KINGDOM
Tel. +44/151.482.4800
Fax +44/151.482.4811

®: Registered Trademark
Published by Terumo Europe N.V.

www.terumo.com
www.terumo-europe.com